


School balance KERN EMB-S

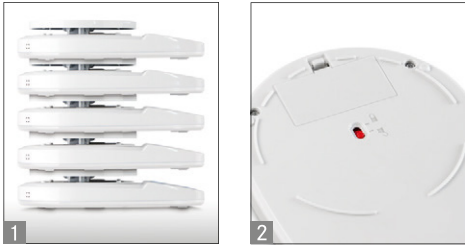

The space miracle – stackable school balance for experiment-based teaching

Features

Technical data

Accessories

- **Simple and convenient 2-key operation**, making them ideal for use in schools and universities
 - **Tare function** facilitates formulation work
 - **Particularly flat design**
 - **Secure and non-slip positioning** with rubber feet
 - **1 Stackable** for space-saving storage
 - **2 A practical transport safety lock** prevents the weighing technology from being overloaded when stored or during transportation
 - **Hook for underfloor weighing** standard
- Large LCD display, digit height 15 mm
 - Dimensions weighing surface Ø 120 mm
 - Weighing plate material plastic
 - Net weight approx. 0,5 kg/piece
 - Optional battery operation, 2 x 1,5 V AA, standard, operating time up to 35 h, AUTO-OFF function to preserve the battery
 - Permissible ambient temperature 5 °C / 35 °C
 - **Note:** The models are only delivered in a set of 5 units. i.e. the price given in the table refers to a delivery of 5 items. Cannot be delivered individually. The calibration prices given here refer to calibration of a single balance

- **External universal mains adapter**, with universal input and optional input socket adapters for EU, GB, USA, CH, not included, can be retrofitted, KERN YKA-03

STANDARD

CAL EXT

UNIT

BATT

DMS

1 DAY

2 YEARS WARRANTY

OPTION

B MULTI

DAkKS +3 DAYS

Model	Weighing range [Max] g	Readout [d] g	Reproducibility g	Linearity g	Dimensions housing WxDxH mm	Option DAkKS Calibr. Certificate	
						DKD KERN	
KERN							
EMB 200-1SS05	200	0,1	0,2	± 0,3	145x205x45	963-127	
EMB 500-1SS05	500	0,1	0,2	± 0,3	145x205x45	963-127	
EMB 2000-OSS05	2000	1	2	± 3	145x205x45	963-127	
EMB 5000-OSS05	5000	1	2	± 3	145x205x45	963-127	

KERN Pictograms:



Internal adjusting: Quick setting up of the balance's accuracy with internal adjusting weight (motordriven).



Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required.



Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



Alibi memory: Electronic archiving of weighing results, complying with the 2009/23/EC standard.



Data interface RS-232: To connect the balance to a printer, PC or network.



RS-485 data interface: To connect the balance to a printer, PC or other peripherals. High tolerance against electromagnetic disturbance.



USB data interface: To connect the balance to a printer, PC or other peripherals.



Bluetooth* data interface: To transfer data from the balance to a printer, PC or other peripherals.



WLAN data interface: To transfer data from the balance to a printer, PC or other peripherals.



Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.



Interface for second balance: For direct connection of a second balance.



Network interface: For connecting the scale to an Ethernet network. With KERN products you can use a universal RS-232/LAN converter.



Wireless data transfer: between the weighing unit and the evaluation unit using an integrated radio module.



GLP/ISO log: The balance displays the weight, date and time, regardless of a printer connection.



GLP/ISO log: With weight, date and time. Only with KERN printers.



Piece counting: Reference quantities selectable. Display can be switched from piece to weight.



Recipe level A: Separate memory for the weight of the tare container and the recipe ingredients (net total).



Recipe level B: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display.



Recipe level C: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, adjustment of recipe when dosages are exceeded, multiplier function, barcode.



Totalising level A: The weights of similar items can be added together and the total can be printed out.



Totalising level C: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, adjustment of recipe when dosages are exceeded, multiplier function, barcode recognition.



Percentage determination: Determining the deviation in % from the target value (100 %).



Weighing units: Can be switched to e.g. non-metric units at the touch of a key. See balance model. Please refer to KERN's website for more details.



Weighing with tolerance range: Upper and lower limiting values can be programmed individually for e.g. dosing, sorting and portioning.



Hold function: (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value.



Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram.



ATEX explosion protection: Suitable for use in hazardous industrial environments, in which there is explosion danger. The ATEX marking is specified for each device.



Stainless steel: The balance is protected against corrosion.



Suspended weighing: Load support with hook on the underside of the balance.



Battery operation: Ready for battery operation. The battery type is specified for each device.



Rechargeable battery pack: Rechargeable set.



Universal mains adapter: with universal input and optional input socket adapters for
A) EU, GB
B) EU, GB, CH, USA
C) EU, GB, CH, USA, AUS



Mains adapter: 230V/50Hz in standard version for EU. On request GB, USA or AUS version available.



Power supply: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request.



Weighing principle: Strain gauge
Electrical resistor on an elastic deforming body.



Weighing principle: Tuning fork
A resonating body is electromagnetically excited, causing it to oscillate.



Weighing principle: Electromagnetic force compensation Coil inside a permanent magnet. For the most accurate weighings.



Weighing principle: Single cell technology
Advanced version of the force compensation principle with the highest level of precision.



Verification possible:
The time required for verification is specified in the pictogram.



DAKkS calibration possible (DKD): The time required for DAKkS calibration is shown in days in the pictogram.



Package shipment: The time required for internal shipping preparations is shown in days in the pictogram.



Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram.



Warranty: The warranty period is shown in the pictogram.

KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2000 kg. In combination with a DAKkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAKkS calibration laboratory today is one of the most modern and best-equipped DAKkS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAKkS calibration of

balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DAKkS calibration of balances with a maximum load of up to 50 t
- DAKkS calibration of weights in the range of 1 mg – 2500 kg
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAKkS calibration certificates in the following languages D, GB, F, I, E, NL, PL

Your KERN specialist dealer: